



Lab No.: C 284/03

Date: 24 July 2003

REPORT

(This report is not to be used for advertising purposes)

on Disinfectant Spray

1 sample received from Magna International Pte Ltd

on 15 July 20 03.

SAMPLE LABEL

One bottle of sample, capped and sealed, bearing the following commercial label,

“ Magna-Legionella-X, 500ml
Viral-Free
Multi-Purpose Disinfectant & Cleaner

Magna-G12, Canada ”

TEST METHOD

Modified Kelsey Sykes Capacity Test Option for Hospital Grade Disinfectant under dirty conditions. Four test organisms, *Escherichia coli* NCTC 8196, *Proteus vulgaris* NCTC 4635, *Pseudomonas aeruginosa* NCTC 6749 and *Staphylococcus aureus* NCTC 4163 were used. The inoculum size of each of the test organisms was not less than 2×10^8 or more than 2×10^9 organisms introduced into the individual test samples of the disinfectant solution. The disinfectant sample was tested neat.

The method is essentially that given by Kelsey & Maurer (Kelsey, J.C and Maurer Isobel, M. Pharmaceutical Journal (UK) 213: 528-230, (1974)). The disinfectant is tested at the use-concentration recommended by the manufacturer. The test consists of challenging the disinfectant with bacterial inoculum, withdrawing a sample after a given time (8 minutes) and culturing the sample in a suitable recovery culture medium. After this sampling, the mixture is again challenged by a second inoculum and after a second interval (18 minutes) is again sampled for culturing. The sample is passed or failed according to the extent of growth shown in the two cultures sampled.

Dr Charles Tang
Sr. Principal Scientific Officer
Pharmaceutical Microbiology Laboratory
Department of Pathology



Lab No.: C 284/03

REPORT

TEST RESULT

SAMPLE: " Magna-Legionella-X, Viral-Free, Multi-Purpose Disinfectant & Cleaner "

TEST CONCENTRATION: Neat

ESCHERICHIA COLI NCTC 8196

<u>INCUBATION TIME</u>	<u>POSITIVE/NEGATIVE CONTROLS</u>	<u>SAMPLING/EXPOSURE TIME</u>		<u>REMARK</u>
		<u>8 MINS</u>	<u>18 MINS</u>	
24 HRS	++ / --	-----	-----	Pass
48 HRS	++ / --	-----	-----	

PROTEUS VULGARIS NCTC 4635

<u>INCUBATION TIME</u>	<u>POSITIVE/NEGATIVE CONTROLS</u>	<u>SAMPLING/EXPOSURE TIME</u>		<u>REMARK</u>
		<u>8 MINS</u>	<u>18 MINS</u>	
24 HRS	++ / --	-----	-----	Pass
48 HRS	++ / --	-----	-----	

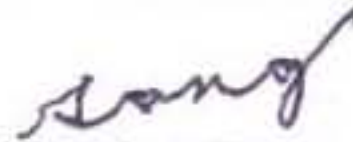
PSEUDOMONAS AERUGINOSA NCTC 6749

<u>INCUBATION TIME</u>	<u>POSITIVE/NEGATIVE CONTROLS</u>	<u>SAMPLING/EXPOSURE TIME</u>		<u>REMARK</u>
		<u>8MINS</u>	<u>18 MINS</u>	
24 HRS	++ / --	-----	-----	Pass
48 HRS	++ / --	-----	-----	

STAPHYLOCOCCUS AUREUS NCTC 4163

<u>INCUBATION TIME</u>	<u>POSITIVE/NEGATIVE CONTROLS</u>	<u>SAMPLING/EXPOSURE TIME</u>		<u>REMARK</u>
		<u>8 MINS</u>	<u>18 MINS</u>	
24 HRS	++ / --	-----	-----	Pass
48 HRS	++ / --	-----	-----	

Note: + : Growth in one tube of recovery broth
- : No growth in one tube of recovery broth


Dr Charles Tang
Sr. Principal Scientific Officer
Pharmaceutical Microbiology Laboratory
Department of Pathology